



## Intersurgical EcoLite™ adult medium concentration oxygen mask



Oxygen and Aerosol Therapy ▪ Oxygen Therapy

## Intersurgical EcoLite™ adult medium concentration oxygen mask

### Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

### Comfortable for the Patient

Improved patient comfort and safety has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC and therefore phthalate free materials in the same product. The polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patients face.

[▶ Watch the video](#)

### Incurved nose seal

Conforms to different nose shapes designed to reduce the risk of oxygen entering patient's eyes

### No metal nose clip

MRI compatible

### Exhalation ports

Designed to minimise the amount of rebreathed carbon dioxide

### A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes


### Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask, reducing the environmental impact.

[Find out more](#)



[Make an enquiry](#)

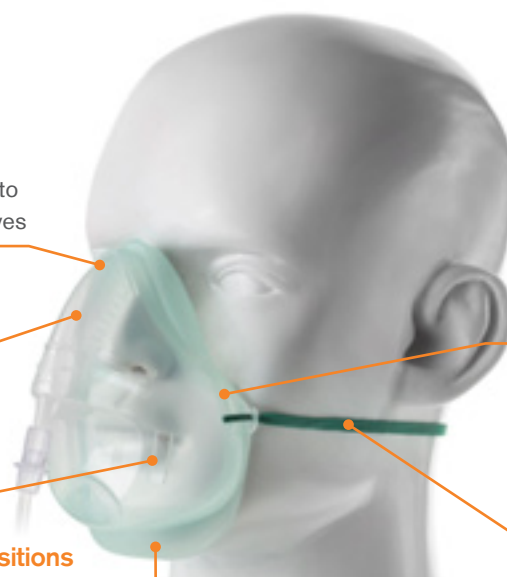

### Lightweight

Mask manufactured using lightweight materials

### Soft face seals

Increased patient comfort

### Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears



Available with oxygen tubing featuring a universal connector, which can be attached to a 6mm flow meter spigot or fitted directly onto a 9/16-18 UNF Diameter Index Safety System (D.I.S.S.) flow meter



### Adult masks

Code	Description	Tube length	Box Qty.
1136015	Intersurgical EcoLite, adult, medium concentration oxygen mask		55
1135015	Intersurgical EcoLite, adult, medium concentration oxygen mask with tube	2.1m	40
1135017	Intersurgical EcoLite, adult, medium concentration oxygen mask and tube with universal connector	2.1m	30

Lower environmental impact product

IS10.16 • Issue 14 02.24



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016, ISO 14001:2015 and MDSAP

**Please think before you print**  
Save energy and paper.  
If you must print this information sheet please print it double sided.